

Instructions for Use

Revolution cfDNA Max 20 Kit™ Revolution cfDNA Max 20 Reagent Kit™

Version E

For in vitro diagnostic use only



REF 402000, 401000



nRichDX, Inc.
15339 Barranca Pkwy
Irvine, CA 92618 USA
+1 949-341-1980

nRich^{DX}
POWERING LIQUID BIOPSY

Table of Contents

Intended Use	3
Summary and Explanation	3
Principles of the Procedure	4
Materials Provided	5
Other Materials Required (Not Provided)	6
Revolution nRicher™ Cartridge Sample Tube Usage and Handling	7
The Revolution Plus Processor and Cartridge Rack	8
Reagent Storing and Handling	9
Sample and Reagent Volumes	9
General Precautions	9
Procedure	
1. Protease Preparation	10
2. Sample Preparation	11
3. nRicher Cartridge Sample Tube Preparation	11
4. Lysis and Bead Binding	12
5. Bead Capture	13
6. Bead Washing	14
7. Elution	16
Troubleshooting	17
Technical Support	18
Warnings and Precautions	19
Appendix	23
Symbols	24
Trademarks, Terms, and Warranty Information	25

Intended Use

The Revolution cfDNA Max 20 Kit, and Revolution cfDNA Max 20 Reagent Kit are intended to extract cell-free DNA (cfDNA) from human plasma and urine.

The kits are intended to be used by physicians and technicians that have received training in molecular biology laboratory techniques.

These kits are intended for in vitro diagnostic use only.

Summary and Explanation

The Revolution cfDNA kits employ well-characterized technology to extract cfDNA from large-volume plasma and urine samples. The kit procedures are designed so users can process multiple samples simultaneously.

The procedures are suitable for nucleic acid isolation from human cell-free plasma or urine. Samples can be either fresh or frozen, although it is recommended that samples that have been previously frozen and thawed are not frozen again. The procedures are designed for minimal user handling, which enables the users to safely handle potentially infectious samples.

The isolated cfDNA is ready for use in downstream applications, including PCR, real-time PCR (RT-PCR), and Next-Generation Sequencing (NGS). Alternatively, the purified cfDNA can be stored at -25 °C to -15 °C for later use.

Principles of the Procedure

Each Revolution cfDNA kit procedure includes the following steps:

- Digest proteins in the plasma or urine sample and protect cfDNA from degradation
- Bind the cfDNA in the plasma or urine sample to magnetic beads
- Capture the beads by magnetic separation
- Wash the beads
- Elute the cfDNA from the magnetic beads

Protease Treatment

Enzymes and other proteins are digested in the Revolution nRicher Cartridge's Sample Tube (see Figure 1 on page 7 for a diagram of the Revolution nRicher Cartridge and its parts).

Addition of Antifoaming Agent

An antifoaming agent is added to each sample to minimize foam formation during the Bead Binding and Bead Capture steps.

Bead Binding

Revolution cfDNA Magnetic Beads and Revolution cfDNA Binding Buffer are combined with the sample in the Revolution nRicher Cartridge's Sample Tube, placed into the Revolution Plus Processor, and incubated to allow the magnetic beads to capture the cfDNA in the sample.

Bead Capture

The Revolution Mag Capsule is attached to the Revolution nRicher Cartridge's Capture Tube, which is again placed into the Revolution Plus Processor and incubated to capture the beads in the Capture Tube portion of the nRicher Cartridge.

Bead Washing

The Capture Tube is removed from the nRicher Cartridge and washed twice with Revolution cfDNA Wash Solution and then rinsed twice with 80% ethanol using the Revolution Mag Rack to capture the magnetic beads after each wash or rinse. The beads are then dried.

Elution

The cfDNA is eluted from the beads with Revolution cfDNA Elution Buffer and is ready for downstream applications.

Materials Provided

IMPORTANT: Upon receipt of the kits, remove the Protease Powder and Magnetic Beads from the kit and store them at the temperatures indicated on the component labels and package insert. All other kit components may be stored at ambient temperature.

Revolution cfDNA Max 20 Kit, (# 402000)

- Revolution cfDNA Surfactant, 7 mL, 414000
- Revolution cfDNA Lysis Buffer, 22 mL, 420000[†]
- Revolution cfDNA Protease Powder, 50 mg, 422000[‡]
- Revolution cfDNA Protease Buffer, 2.8 mL, 418000
- Revolution cfDNA Binding Buffer, 160 mL, 410000[†]
- Revolution cfDNA Antifoaming Agent, 500 μ L, 426000
- Revolution cfDNA Magnetic Beads, 2 x 1.4 mL, 424000[‡]
- Revolution cfDNA Wash Solution, 54 mL, 412000
- Revolution cfDNA Elution Buffer, 3.5 mL, 416000
- Revolution nRicher Cartridge Sample Tubes, 3 packs of 8 Sample Tubes, 450007
- Revolution nRicher Cartridge Capture Tubes, 3 packs of 8 Capture Tubes, 450003
- Revolution Capture Tube Caps, Qty. 100, 400063
- Instructions for Use (This document. The most current version is available at www.nrichdx.com)

Revolution cfDNA Max 20 Reagent Kit, (# 401000)

- Revolution cfDNA Surfactant, 7 mL, 414000
- Revolution cfDNA Lysis Buffer, 22 mL, 420000[†]
- Revolution cfDNA Protease Powder, 50 mg, 422000[‡]
- Revolution cfDNA Protease Buffer, 2.8 mL, 418000
- Revolution cfDNA Binding Buffer, 160 mL, 410000[†]
- Revolution cfDNA Antifoaming Agent, 500 μ L, 426000
- Revolution cfDNA Magnetic Beads, 2 x 1.4 mL, 424000[‡]
- Revolution cfDNA Wash Solution, 54 mL, 412000
- Revolution cfDNA Elution Buffer, 3.5 mL, 416000
- Instructions for Use (This document. The most current version is available at www.nrichdx.com)

[†] Contains guanidinium thiocyanate. Do not combine with disinfectants that contain bleach. See “Warnings and Precautions” section of this document for more information.

[‡] Upon receipt, store the Protease Powder frozen at -25 °C to -15 °C and the Magnetic Beads at the temperature indicated on the Magnetic Beads label and Package Insert. All other kit components may be stored at ambient temperature 15° to 25° C (59° to 77° F).

Other Materials Required (Not Provided)

Always wear personal protective equipment, such as a lab coat, protective eyewear, and disposable gloves when working with chemicals. Consult the appropriate Safety Data Sheets (SDSs for Revolution kits are available at <https://www.nrichdx.com/material-data-safety-sheets>; other SDSs are available from the product supplier) for more information on safe handling and use.

Revolution Sample Prep System Equipment

- Drip Tray, 100291
- Revolution Cartridge Rack, 200600
- Revolution Mag Capsules, 200700
- Revolution Capture Tube Mag Rack, 200800
- Revolution Plus Processor, 2000-PLUS

Additional Materials and Reagents

- 80% ethanol
- Centrifuge and Microcentrifuge[†]
- Non-magnetic microvial rack* (**NOTE:** Capture Tube will generally fit a microvial tube rack)
- Phosphate Buffer Saline (PBS) pH 7.5
- Pipettors[†], pipet tips[‡], and serological pipettes
- Vortex instrument with 2mL microvial tube adaptor*

***NOTE:** The nRicher Cartridge's Capture Tube's width will typically fit into a standard non-magnetic microvial rack or vortex instrument microvial tube adaptor.

nRicher Cartridge Usage and Handling

The nRicher Cartridge is comprised of two parts - the Sample Tube and the Capture Tube as shown in Figure 1. To add the Capture Tube push the Capture Tube's open end evenly into the open end of the Sample Tube until you hear a click sound and physical sensation as shown in Figure 2. To remove the Capture Tube pull the Capture Tube evenly away from the Sample Tube as shown in Figure 2.

[†] We strongly recommend that instruments are calibrated at regular intervals to ensure that samples are processed consistently and accurately

[‡] We strongly recommend using pipette tips with aerosol barriers to prevent cross contamination.

Similarly, the Revolution Mag Capsule is added to the Capture Tube by placing the Mag Capsule's larger opening over the top of the Capture Tube and pushing down until the Mag Capsule is fully seated on the Capture Tube as shown in Figure 3. Note: a slight clockwise/counter-clockwise twisting of the Mag Capsule until the Mag Capsule is fully seated. To remove the Mag Capsule pull the Mag Capsule evenly away from the Capture Tube as shown in Fig 3.

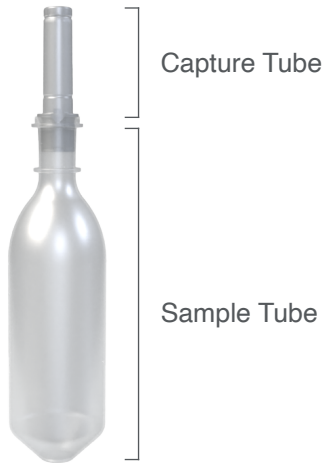


Figure 1 - nRicher Cartridge
The nRicher Cartridge is comprised of the Sample Tube and the Capture Tube

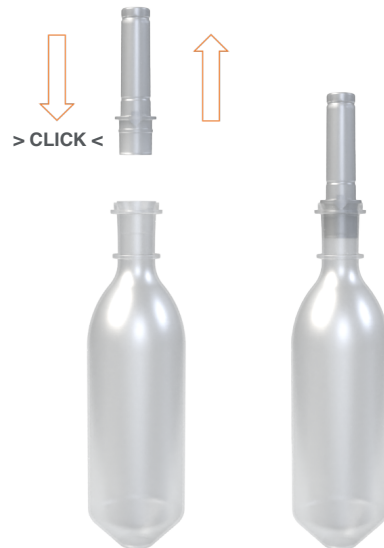


Figure 2 - Capture Tube
Joining or separating the Capture Tube and Sample Tube

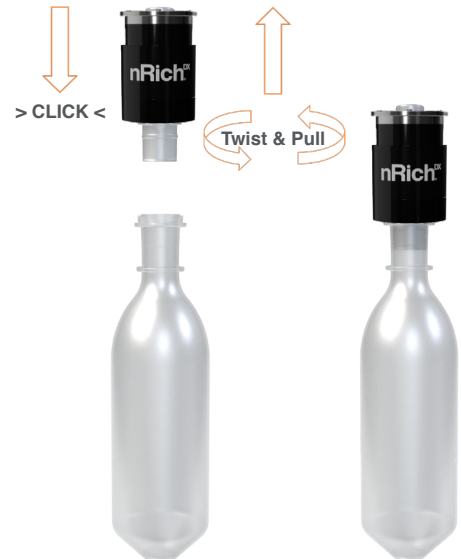


Figure 3 - Mag Capsule
Joining or separating the Mag Capsule and Capture Tube from the Sample Tube



Figure 4 - Releasing the Capture Tube from the Mag Capsule
Place the Mag Capsule with the open end of the Capture Tube facing upward on a laboratory benchtop and press down; the Capture Tube will click as it moves upward slightly and is released from the Mag Capsule. The Capture Tube may now be removed from the Mag Capsule by gripping the Capture Tube and pulling it upward gently away from the Mag Capsule.



Figure 5 - Magnetic Rack
12-position Capture Tube Magnetic Rack (Mag Rack)



Figure 6 - Cartridge Rack with Cartridges
nRicher Cartridges shown in the Cartridge Rack; the Capture Tube is removed and placed in a separate rack when accessing the Sample Tube portion of the nRicher Cartridge.



Figure 7 - Rack with Cartridges, Mag Capsules
nRicher Cartridges with Mag Capsules attached and in the Cartridge Rack.



Figure 8
The Revolution Plus Processor.



Figure 9
The Revolution Plus Processor - door open. Note the Cartridge Rack attachment rod and pins - two pins on the left side of the rod and one pin on the right side. Pins attach the Cartridge Rack to the Processor



Figure 10 - Placing the Cartridge Rack into and removal from, the Revolution Plus Processor

Two pin slots on the left-hand side of the Cartridge Rack and one pin slot on the right-hand side securely attach the rack to the processor as follows: Grip the Cartridge Rack handles and place the Cartridge Rack directly over the attachment rod so the pins are to the left and adjacent to the slots. Slide the rack to the left into the pins on the attachment rod until the pins are attached to the Cartridge Rack and the single pin is positioned at the red arrow shown above. To remove the Cartridge Rack, first ensure all motion of the Rack has stopped, grip the Cartridge Rack handles and slide the rack to the right until the pins have detached from the Cartridge Rack's slots. Lift the rack over the attachment rod and free from the Revolution Plus Processor.

Reagent Storing and Handling

The Revolution cfDNA kits are shipped at room temperature.

IMPORTANT: Upon arrival, remove the components indicated below and store them at the indicated storage temperatures:

- Revolution cfDNA Protease Powder should be stored at -25 °C to -15 °C and can be used until the kit expiration date without affecting protease performance.
- Revolution cfDNA Magnetic Beads should be stored at the temperature indicated on the label and the package insert.
- All other kit components can be stored at ambient temperature (15 °C to 30 °C) until the expiration date without affecting component performance.

Sample and Reagent Volumes

The reagent volumes in the protocol differ based on the starting sample volume. Tables with all of the sample and reagent volumes are included in the Appendix. To make the protocol easy to use, the following icons will be used throughout the protocol to indicate the correct reagent volume for each starting sample volume.

- Green square (■) for **5 mL** sample volume
- Blue circle (●) for **10 mL** sample volume
- Tan diamond (◆) for **15 mL** sample volume
- Red triangle (▲) for **20 mL** sample volume

For example, ■ 250 μ L indicates that at this particular step, 250 μ L would be added if your starting volume was 5 mL.

General Precautions

- Perform all steps at ambient temperature (15 °C to 30 °C) unless otherwise noted.
- If you observe a precipitate in the Revolution cfDNA Lysis Buffer, incubate the Lysis Buffer

- at 37 °C until the precipitate dissolves. This can occur if storage temperatures are too low.
- If you observe a precipitate in the Revolution cfDNA Surfactant, incubate the Surfactant at 37 °C until the precipitate dissolves. This can happen if storage temperatures are too low.
 - Before use, thoroughly mix the Revolution Antifoaming Agent by vortexing at high speed for 10 seconds. As it is a suspension, complete mixing is essential for consistent performance. Due to its high viscosity, exercise caution and dispense slowly when pipetting the Antifoaming Agent.
 - Vortex the Revolution cfDNA Magnetic Beads to fully resuspend them immediately before use.
 - Blood samples should be collected in K2EDTA tubes or Streck Cell-Free DNA BCT tubes.
 - Cell-free DNA is stable for 24 hours in K2EDTA plasma tubes stored at 2 °C to 10 °C.
 - Blood samples collected in Streck Cell-Free DNA BCT tubes remain stable for up to 14 days due to the formaldehyde-free preservative in the tubes.
 - For urine samples, use fresh, preferably first-void urine, stored at 2 °C to 10 °C and spun down to cell-free status within 24 hours of collection to avoid increases in genomic DNA (gDNA) and microbial growth.
 - If it is not possible to process urine samples immediately after collection, cfDNA urine preservative should be added. Samples collected in the preservative should preferentially be stored and centrifuged per manufacturer recommendations.
 - Adequate Capture Tube caps are included in the kit to ensure a new cap is used for each capping step in the protocol. To prevent cross-contamination do not reuse caps.

Procedure

1. Protease Preparation

- 1.1. Transfer 2.5 mL of Protease Buffer (PB) to the bottle of Protease Powder (PP).
- 1.2. Cap the bottle and gently invert 8 to 10 times to dissolve the powder.
- 1.3. Place the rehydrated protease solution on ice until use.

NOTE: Unused rehydrated protease can be stored at 2 °C to 10 °C until needed for re-use.

2. Sample Preparation

For Plasma Samples

- 2.1. Centrifuge the blood samples at 2,000 x g for 10 minutes at 2 °C to 10 °C.
- 2.2. Transfer the plasma to a new centrifuge tube.
- 2.3. Centrifuge the plasma samples at 16,000 x g for 10 minutes at 2 °C to 10 °C.
NOTE: Alternatively, the plasma can be centrifuged at 6,000 x g for 30 minutes to remove any residual blood and cell debris.
- 2.4. Transfer the cell-free supernatants into fresh tubes.
- 2.5. All samples must be equilibrated to one of the starting sample volumes (e.g., 5, 10, 15, 20 mL) based on the following criteria:
 - Samples < 5 mL — bring the volume to ■ 5 mL with PBS pH 7.5
 - Samples < 10 mL — bring the volume to ● 10 mL with PBS pH 7.5
 - Samples < 15 mL — bring the volume to ◆ 15 mL with PBS pH 7.5
 - Samples < 20 mL — bring the volume to ▲ 20 mL with PBS pH 7.5
- 2.6. Process samples immediately or store on ice until use.

For Urine Samples

- 2.7. Centrifuge the urine samples at 16,000 x g for 10 minutes at 2 °C to 10 °C.
- 2.8. Transfer the cell-free supernatants into fresh tubes.
- 2.9. All samples must be equilibrated to one of the starting sample volumes (e.g., 5, 10, 15, 20 mL) based on the following criteria:
 - Samples < 5 mL — bring the volume to ■ 5 mL with PBS pH 7.5
 - Samples < 10 mL — bring the volume to ● 10 mL with PBS pH 7.5
 - Samples < 15 mL — bring the volume to ◆ 15 mL with PBS pH 7.5
 - Samples < 20 mL — bring the volume to ▲ 20 mL with PBS pH 7.5
- 2.10. Process samples immediately or store on ice until use.

3. nRicher Cartridge Sample Tube Preparation

- 3.1. Label the nRicher Cartridge elements (both the Sample Tube and Capture

Tube; see Figure 1, page 7) with a sample identifier; use one Sample Tube per sample.

- 3.2. Place the Sample Tube(s) into the Cartridge Rack

NOTE: The Sample Tubes should be left in the rack for the entirety of the extraction.

4. Lysis and Bead Binding

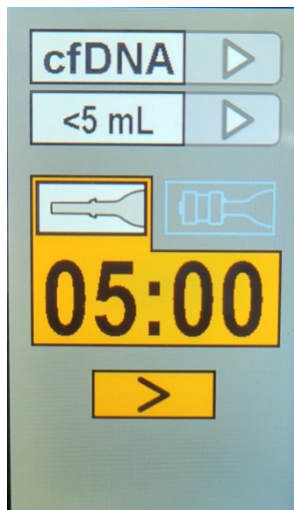




Figure 11 - GUI for bead binding on the Revolution Plus

- 4.1. Add cell-free plasma/urine to the labeled Sample Tube.
- 4.2. Add ■ 250 μL , ● 500 μL , ◆ 750 μL , or ▲ 1000 μL Surfactant Solution to each Sample Tube.
- 4.3. Add ■ 800 μL , ● 1600 μL , ◆ 2400 μL , or ▲ 3200 μL Lysis Buffer to each Sample Tube.
- 4.4. Add ■ 80 μL , ● 160 μL , ◆ 240 μL , or ▲ 320 μL Protease Solution to each Sample Tube.
- 4.5. Add ■ 6 mL, ● 12 mL, ◆ 18 mL, or ▲ 24 mL Binding Buffer to each Sample Tube.
- 4.6. Vortex the Anti-Foaming Agent for 20 seconds at high speed. For all sample volumes, add 10 μL Antifoaming Agent to each sample tube.
- 4.7. Resuspend the Revolution cfDNA Magnetic Beads by vortexing at medium


speed for 10 seconds.

- 4.8. Add ■ 100 μL , ● 100 μL , ◆ 150 μL , or ▲ 150 μL Magnetic Beads to each Sample Tube.
- 4.9. Attach the Capture Tube to its companion Sample Tube, and then transfer the Cartridge Rack containing the assembled nRicher Cartridges to the Revolution Processor using the holes positioned on the rack (see Figures 9 and 10, on page 8).
- 4.10. Choose the  icon under “cfDNA” and adjust the sample volume accordingly on the Revolution Plus’ Graphical User Interface (GUI). Initiate the process by pressing the  icon.
- 4.11. When the Revolution Plus Processor stops and illuminates orange, remove the Cartridge Rack from the device (see Figure 10, page 8), and place the rack on a level surface.

5. Bead Capture



Figure 12 - GUI for bead bead capture on the Revolution Plus

- 5.1. Snap a Mag Capsule onto the Capture Tube of each nRicher Cartridge (Figure 3, page 7), and then place the Cartridge Rack into the Revolution Plus Processor using the holes positioned on the rack to securely attach the Cartridge Rack to the Revolution Plus Processor (Figures 9 and 10, page 8).
- 5.2. Choose the  icon under “cfDNA” and adjust the sample volume accordingly on the Revolution Plus’ GUI. Initiate the procedure by pressing

the  icon.

- 5.3. When the Revolution Plus's motion stops and the GUI illuminates blue, remove the Cartridge Rack from the processor, and place the rack on a level surface; let the rack sit for 1 minute to allow all of the liquid to drain from the nRicher Cartridge's Capture Tube into the Sample Tube.
- 5.4. Remove the Mag Capsule and Capture Tube from each Sample Tube by twisting the Mag Capsule counterclockwise (Figure 3, page 7).
IMPORTANT: Ensure the Mag Capsule remains attached to the Capture Tube for this step.
- 5.5. Gently orient the Mag Capsule and Capture Tube so that the open end of the Capture Tube is facing upward; then press the Mag Capsule down toward the benchtop to release the Capture Tube from the Mag Capsule as shown in Figure 4, page 7. Place the released Capture Tube(s) in a non-magnetic microvial rack.
- 5.6. Discard the liquid remaining in each Sample Tube along with the Sample Tube itself in biohazardous waste.

6. Bead Washing

- 6.1. Add 1 mL Wash Solution to each Capture Tube.
- 6.2. Firmly seal each Capture Tube with a new cap and vortex for 10 seconds at medium speed. **IMPORTANT:** Ensure the cap is pressed firmly into the capture tube prior to vortexing to prevent leakage.
- 6.3. Inspect Capture Tube(s) to ensure beads are fully resuspended; if not, then vortex for an additional 10 seconds or until the beads are fully resuspended (no bead clumping is observed).
- 6.4. Centrifuge the Capture Tube(s) for 1 second in a microcentrifuge to collect all remaining liquid from the cap to the bottom of the tube.
- 6.5. Place the Capture Tube(s) into the Mag Rack (Figure 5, page 7) for 2 minutes.
- 6.6. Inspect the Capture Tube(s) to ensure that the supernatant is clear before proceeding to the next step; if the supernatant is not clear, leave the

- Capture Tube(s) on the Mag Rack for an additional 1 minute.
- 6.7. Remove the cap from each Capture Tube and carefully aspirate and dispose of the supernatant from each tube while the Capture Tube(s) remain in the Mag Rack. **CAUTION:** Be careful not to disturb the bead line and to use a sterile pipette tip for each sample.
 - 6.8. Transfer the Capture Tube(s) from the Mag Rack to a non-magnetic microvial rack.
 - 6.9. Repeat steps 6.1 to 6.7 for a second wash with Wash Solution.
 - 6.10. Add 1 mL 80% ethanol rinse to each Capture Tube.
 - 6.11. Recap the Capture Tube(s) and vortex at medium speed for 10 seconds.
 - 6.12. Centrifuge the Capture Tube(s) for 1 second in a microcentrifuge to collect all remaining liquid from the cap to the bottom of the Capture Tube.
 - 6.13. Place the Capture Tube(s) in the Mag Rack for 2 minutes.
 - 6.14. Inspect the Capture Tube(s) to ensure that the supernatant is clear before proceeding to the next step; if the supernatant is not clear, leave the Capture Tube(s) on the Mag Rack for 1 minute.
 - 6.15. Remove the cap from each Capture Tube. Carefully aspirate and dispose of the supernatant from each Capture Tube (using a sterile pipette tip for each sample) while the Capture Tube(s) remain in the Mag Rack.
 - 6.16. Transfer the Capture Tube(s) to a non-magnetic microvial rack. **NOTE: do not repeat this step after the second ethanol rinse. Go to step 6.18.**
 - 6.17. Repeat rinse steps 6.10 to 6.15 for a second ethanol rinse.
 - 6.18. After step 6.15, tap the Mag Rack on the bench 5 times to collect all remaining ethanol to the bottom of the Capture Tube. Carefully aspirate any remaining ethanol from the walls and bottom of the Capture Tube to ensure that the drying time in step 6.20 effectively removes any remaining ethanol.
 - 6.19. For all sample volumes, pulse centrifuge the Capture Tube(s) for 1 second in a microcentrifuge to collect all the beads to the bottom of the Capture Tube.
 - 6.20. **BEAD DRYING:** Dry the Capture Tube(s) in a 45 °C heat block for 20 minutes. After drying is complete, transfer the Capture Tube(s) to a non-magnetic microvial tube rack.

NOTE: It is important that the magnetic bead pellet is completely dry and

free of ethanol before proceeding. If the beads are not dry enough, place the capture tubes in the 45°C heat block for 1 to 2 more minutes; however, do not overdry the magnetic pellet which may cause the pellet to be difficult to resuspend and lower overall cfDNA yield.

7. Elution

- 7.1. Add 50 μ L Elution Buffer to each Capture Tube, and recap.
NOTE: Elution buffer volume may range from 25 μ L to 100 μ L. The cfDNA concentration will decrease accordingly when eluting in larger volumes.
- 7.2. Briefly vortex the Capture Tube(s) at medium to high speed for 5 to 10 seconds.
- 7.3. Make sure the magnetic beads are completely resuspended in the elution buffer; if the beads are not completely resuspended (or if bead clumping is observed), vortex for another 5 to 10 seconds.
- 7.4. Vortex the Capture Tube(s) on a vortex shaker for 5 minutes at medium to high speed.
- 7.5. Centrifuge the Capture Tube(s) for 1 second in a microcentrifuge to collect all remaining liquid from the cap to the bottom of the tube.
- 7.6. Place the Capture Tube(s) into the Mag Rack for 2 minutes.
- 7.7. Inspect the Capture Tube(s) to ensure that the supernatant is clear before proceeding to the next step; if the supernatant is not clear, leave the Capture Tube(s) on the Mag Rack for another 2 minutes.
- 7.8. Transfer the eluate into a clean, properly labeled tube or plate, making sure not to disrupt the magnetic bead pellet while collecting the eluate (*you may transfer the eluate to any preferred tube or plate suitable for nucleic acid, the elution volume, and your downstream method; for example, a microvial tube, PCR tube, or 96-well plate*).
- 7.9. Discard the Capture Tube(s) containing magnetic beads as biohazardous waste.
- 7.10. Store the eluate at -25 °C to -15 °C until ready for downstream analysis.

[End of the nRichDX Revolution Plus cfDNA extraction protocol]

Troubleshooting

Observation	Possible cause	Recommended action	
Lower yield than expected	The Revolution cfDNA Magnetic Beads were not properly stored	Store the Revolution cfDNA Magnetic Beads at the temperature indicated on the label and package insert. Do not freeze the beads.	
	An insufficient amount of Revolution cfDNA Magnetic Beads was added	Immediately before use, vortex the tube containing the magnetic beads thoroughly until fully resuspended.	
	The Revolution cfDNA Magnetic Beads are not optimally dried		Drying times may vary depending on the amount of beads used and the environment. Lower volumes of beads require less time to dry. Airflow and humidity in the immediate environment may shorten or lengthen the optimal bead drying time.
			Overdried beads will stick to the wall of the tube and be difficult to resuspend. Gently scrape the beads off the plastic wall using a pipette tip.
			Underdried beads may carry ethanol into the eluate and negatively impact downstream applications. Dry beads slightly longer (1-minute intervals) and make note of the optimal drying time for the specific volume.
	Nucleases	Nuclease contamination will result in a lower yield of intact cfDNA and/or cfRNA. Ensure reagents, pipet tips and plasticware in direct contact with the sample are free of undesired nucleases. Use nuclease free barrier filter tips and filtered hood to minimize presence of airborne nucleases.	
	Magnetic bead clumping is observed	Vortex the tube containing the Magnetic Beads until they are fully resuspended	
The sample contains low levels of cfDNA	Increase the starting sample volume.		
Magnetic bead carryover	Loose beads present in the eluate or inadvertently transferred	Be sure to leave the Capture Tube(s) on the Mag Rack when removing the eluate containing the cfDNA.	
		If beads are carried over into the new tube, place the tube on the Mag Rack again, wait for the beads to pellet and then transfer the sample to another tube.	

Observation	Possible cause	Recommended action
Fluid leaking - Capture Tube not vertical and/or not snug on Sample Tube	Capture Tube not connected correctly to sample tube	Unscrew Capture Tube from Sample Tube. Reattach the Capture Tube to the Sample Tube using the following steps. Lower the Capture Tube opening over the Sample Tube opening. Apply medium pressure to the connection point until a click is heard or sensing the Capture Tube is attached. If the Capture Tube doesn't click, then while applying this pressure, rotate the Capture Tube slightly clockwise and counter-clockwise until the click is heard or sense the tube is attached.
Variations in cfDNA yield from donor to donor	Variation in amount of circulating cfDNA. Levels of cfDNA in circulation can range from 1 to 100 ng/mL of plasma or serum depending on the donor.	For samples containing low levels of cfDNA, increase the starting sample volume.

Technical Support

For additional questions, please contact technical support services at technicalsupport@nrchdx.com

Warnings and Precautions

For In Vitro Diagnostic Use

Users should wear personal protective equipment as required by local laboratory procedures when performing an isolation, including a lab coat, protective eyewear, and disposable nitrile gloves (or equivalent). Please refer to the relevant safety data sheets (SDSs) for more information.

Discard all used materials as biohazardous waste according to local regulations.

CAUTION: The Revolution Lysis Buffer contains guanidinium thiocyanate, which when combined with bleach, forms highly reactive compounds.



CAUTION: DO NOT directly add bleach or acidic solutions to the isolation waste.

Clean up all spills with appropriate laboratory-grade detergent and water. Any spills that contain potentially infectious materials should be cleaned first with laboratory detergent and water followed by 1% (v/v) sodium hypochlorite.

If any of the reagent bottles or containers are damaged and leaking fluids, wear gloves and protective eyewear when discarding the bottles.

Revolution cfDNA Surfactant

Hazard pictograms:




Signal word: DANGER

Hazard and precautionary statements:

H315: Causes skin irritation; H318: Causes serious eye damage; P264: Wear protective gloves / protective clothing / eye protection / face protection; P332 + P313: If skin irritation occurs: Get medical advice/attention; P302 + P352: IF ON SKIN - Wash with plenty of soap and water; P305 + P351 + P338: IF IN EYES - Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing; P333 + P313: If skin irritation or rash occurs, get medical advice/attention; P310: Immediately call a POISON CENTER or doctor / physician.

Revolution cfDNA Lysis Buffer

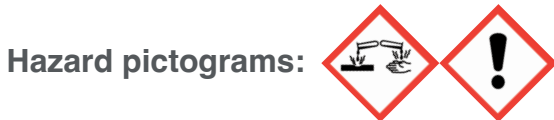

Signal word: DANGER

Hazard and precautionary statements:

H302 + H312 + H332: Harmful if swallowed, in contact with skin or inhaled; H314: Causes severe skin burns and eye damage; H412: Harmful to aquatic life with long lasting effects; P261: Avoid breathing dust / fumes / gas / mist / vapors / spray; P264: Wash skin thoroughly after handling; P270: Do not eat, drink or smoke when using this product; P271: Use only outdoors or in well-ventilated area; P273: Avoid release to the environment; P280: Wear protective gloves / protective clothing / eye protection / face protection; P301 + P312 + P330: IF SWALLOWED - Call a POISON CENTER / doctor if you feel unwell. Rinse mouth; P301 + P330 + P331: IF SWALLOWED - Rinse mouth. Do NOT induce vomiting; P303 + P361 + P353: IF ON SKIN (or hair): Take off immediately all contaminated clothing. rinse skin with water / shower; P304 + P340 + P310: IF INHALED - Remove person to fresh air and keep comfortable for breathing. Immediately call a POISON CENTER / doctor; P305 + P351 + P338 + P310: IF IN EYES - Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing. Immediately call a POISON CENTER / doctor; P363: Wash contaminated clothing before reuse; P405: Store locked up; P501: Dispose of contents / container in an approved waste disposal plant.

Revolution cfDNA Protease Powder




Signal word: DANGER

Hazard and precautionary statements:

H315: Causes skin irritation; H319: Causes serious eye irritation; H334: May cause allergy or asthma symptoms or breathing difficulties if inhaled; H335: May cause respiratory irritation; P261: Avoid breathing dust / fumes / gas / mist / vapors / spray; P305 + P351 + P338: IF IN EYES - Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing; P341 + P311: If experiencing respiratory symptoms: Call a POISON CENTER or doctor / physician.

Revolution cfDNA Protease Buffer

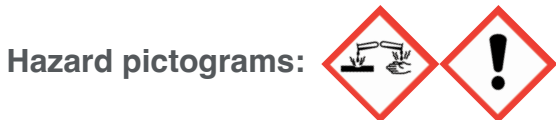

Signal word: DANGER

Hazard and precautionary statements:

H315: Causes skin irritation; H318: Causes serious eye damage; P280: Wear protective gloves / protective clothing / eye protection / face protection; P264: Wash hands thoroughly after handling; P332 + P313: If skin irritation occurs, get medical advice / attention; P302 + P352: IF ON SKIN - Wash with plenty of soap and water; P305 + P351 + P338: IF IN EYES - Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing; P333 + P313: If skin irritation or rash occurs, get medical advice / attention; P310: Immediately call a POISON CENTER or doctor / physician.

Revolution cfDNA Binding Buffer



Signal word: DANGER

Hazard and precautionary statements:

H302: Harmful if swallowed; H314: Causes severe skin burns and eye damage; H412: Harmful to aquatic life with long lasting effects; P260: Do not breathe dust / fume / gas / vapors / spray; P264: Wash hands thoroughly after handling; P270: Do not eat, drink or smoke when using this product; P273: Avoid release to the environment; P280: Wear protective gloves / protective clothing / eye protection / face protection; P301 + P310: IF SWALLOWED - Immediately call a POISON CENTER or doctor / physician; P303 + P361 + P353: IF ON SKIN (or hair) - Remove / take off immediately all contaminated clothing. Rinse skin with water / shower; P 305 + P351 + P338: IF IN EYES - Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing; P310: Immediately call a POISON CENTER or doctor / physician; P330: Rinse mouth.

Revolution cfDNA Wash Solution

Hazard pictograms:



Signal word: WARNING

Hazard and precautionary statements:

H302: Harmful if swallowed; H315: Causes skin irritation; H319: Causes serious eye irritation; P264: Wash hands thoroughly after handling; P270: Do not eat, drink or smoke when using this product; P280: Wear protective gloves / protective clothing / eye protection / face protection; P330: Rinse mouth; P332 + P313: If skin irritation occurs, get medical advice / attention; P305 + P351 + P338: IF IN EYES: Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing. P301 + P310: IF SWALLOWED: Immediately call a POISON CENTER or doctor / physician; P302 + P352: IF ON SKIN: Wash with plenty of soap and water.

Appendix

Table 1. Sample and Reagent Volumes

Sample volume (mL)	Surfactant (μL)	Lysis buffer (μL)	Protease solution (μL)	Binding buffer (mL)	Antifoaming Agent (μL)	Magnetic beads (μL)	Total volume (mL)
5	250	800	80	6	10	100	12.2
10	500	1600	160	12	10	100	24.4
15	750	2400	240	18	10	150	36.6
20	1000	3200	320	24	10	150	48.7

Table 2. Bead Binding and Bead Capture Incubation Times

Sample volume (mL)	Binding Incubation Time (minutes)	Bead Capture Incubation Time (minutes)
5	5	20
10	15	20
15	15	45
20	15	45

*Total volume does not include wash solution

Final elution volume is in a range of 25 μL - 100 μL ; 50 μL is the default recommended elution volume

Symbols



In vitro diagnostic medical device



Catalog numbers



Manufacturer



Use-by date



Batch code



Consult instructions for use



Caution



Temperature range

Trademarks, Terms, and Warranty Information

Trademarks: including the nRichDX® logo and taglines, nRicher™, Revolution Plus™, Revolution Pro™, Revolution Sample Prep System™, Revolution cfTNA Max 20 Kit™, Revolution cfTNA Max 20XL Kit™, Revolution cfTNA Max 50 Kit™, Revolution cfTNA Reagent Kit™, Revolution cfTNA Max 20 Reagent Kit™, Revolution cfTNA Max 50 Reagent Kit™, Revolution cfDNA Max 20 Kit™, Revolution cfDNA Max 20XL Kit™, Revolution cfDNA Max 50 Kit™, Revolution cfDNA Reagent Kit™, Revolution cfDNA Max 20 Reagent Kit™, Revolution cfDNA Max 50 Reagent Kit™, “Powering Liquid Biopsy”, “Driving Diagnostic Excellence®,” and “Liquid Biopsy is revolutionizing cancer treatment. nRichDX empowers that revolution.®,” are trademarks of nRichDX, Inc.

nRichDX products are covered by one or more of the following patents: 10,329,554; 10,585,101; 10,927,366; D886,320; D890,362; Canada patent numbers 3,016,895; 3,042,426; 190733; Europe patent numbers 001720118; 3535058; China patent number ZL2017800688437; Japan patent number 389863; South Korea patent numbers 30-1062125; 10-2239349; Australia patent numbers 201916131; 2017355635, with other patents pending.

Use of this product indicates the agreement of purchaser or user of this kit to the following terms:

1. This product can only be used in accordance with the protocols and procedures provided in this document using the components in this kit. nRichDX, Inc. grants no license under its intellectual property to employ or integrate the components of their kits with any components not included in the kit except as indicated in the protocols and procedures described herein.
2. nRichDX, Inc. does not warranty that this kit or the use of this kit does not infringe upon the rights of third-party entities.
3. This kit and the components therein are only licensed for one-time use; the kit or components cannot be reused, refurbished or resold.
4. nRichDX disclaims any licenses, expressed or implied, other than the licenses that are expressly stated.
5. The purchaser and the user of this kit agree to not act or permit anyone else to act in any way that might lead to or enable any of the prohibited acts indicated above. nRichDX may enforce the prohibitions of this Limited License Agreement in Court and shall work to recover all investigative and Court costs, including attorney fees, in any action that nRichDX takes to enforce this Limited License Agreement or any intellectual property rights that relate to this kit and/or components of the kit.

© 2024 nRichDX, Inc.

The nRichDX Revolution System™

Liquid Biopsy is revolutionizing cancer treatment.
nRichDX empowers that revolution.®

